

RuiduLance™—More than 10 years of guardianship

Product Features

Safety

Gamma irradiation ensures 5 years of shelf life / Needle tip is sealed in the protective cap / Needle tip automatically retracts after use

Low pain level



US&CN patented low pain level design

Convenient use

Easy to use and carry

Various gauges

21G, 23G, 26G, 28G, 30G, various gauges, meet different kind of skin type

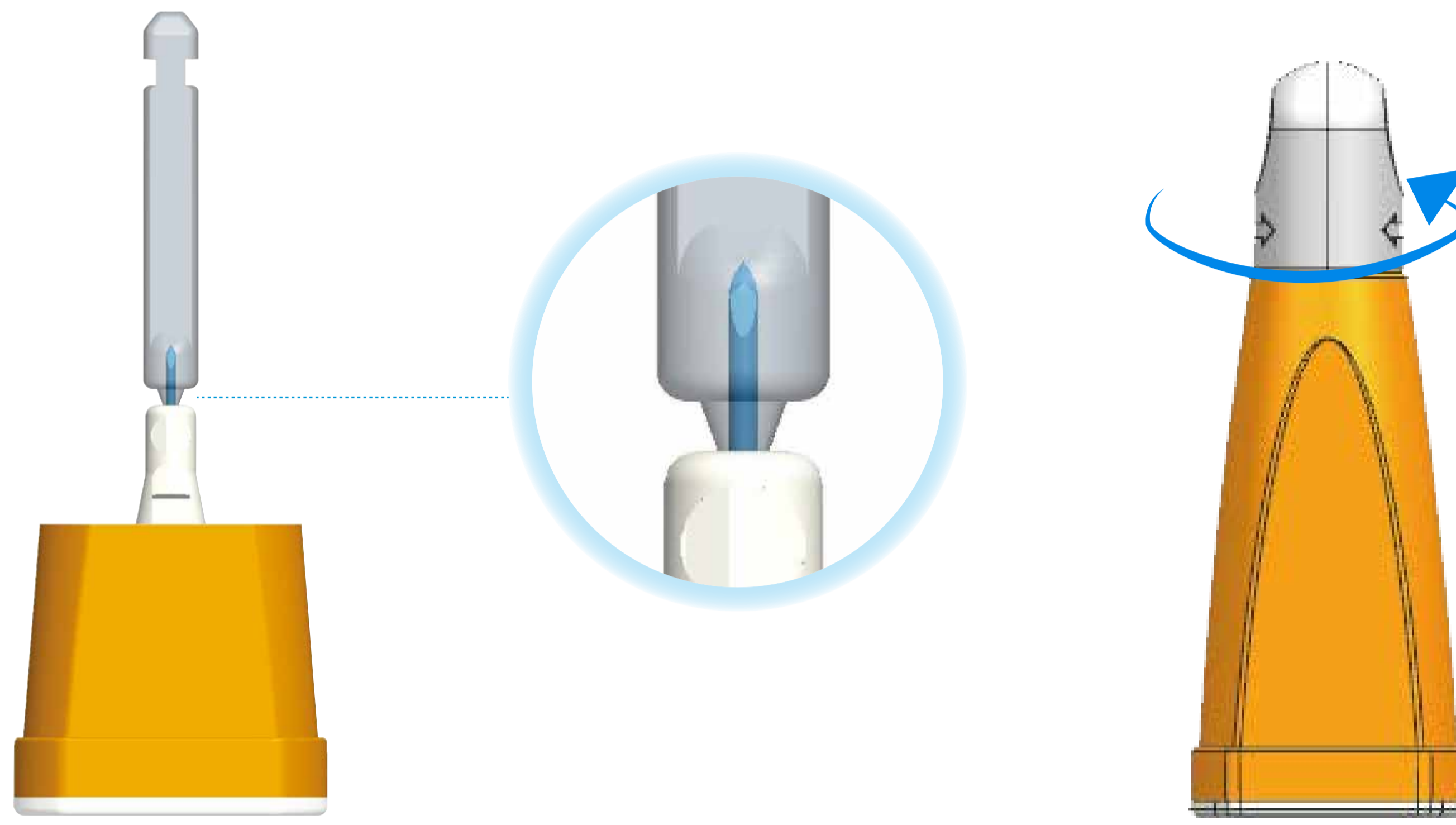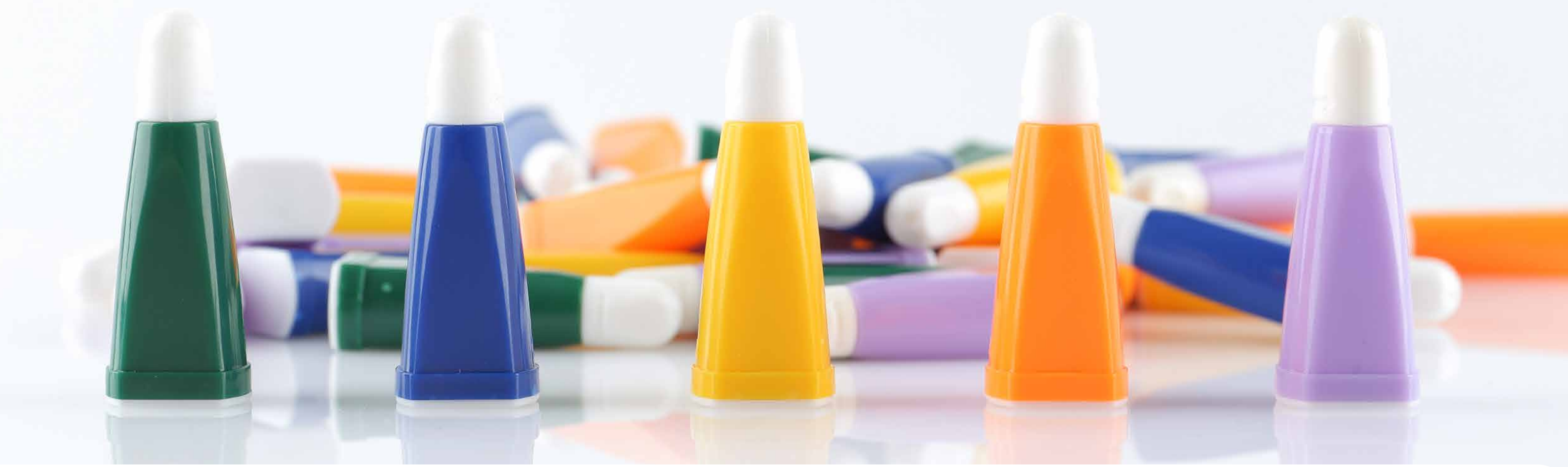

Needle tip is sealed in the plastic handle, ensuring the sterility
Twist off the protective cap, do not pull



US&CN patented design, stable firing system



Invention Patent


US: 9.826.927

CN: ZL 2012 1 0080546.3


- Stable firing system
- Elastic inner sleeve
- Actual-location puncture
- No repeated puncture

Code	Specification	Color	Blood	Package
I17A2122100	21G*2.2mm			100 pcs/box 20 boxes/carton (or per customers' requests)
I17A2322100	23G*2.2mm			
I17A2618100	26G*1.8mm			
I17A2818100	28G*1.8mm			
I17A3016100	30G*1.6mm			


Use Instructions for RuiduLance

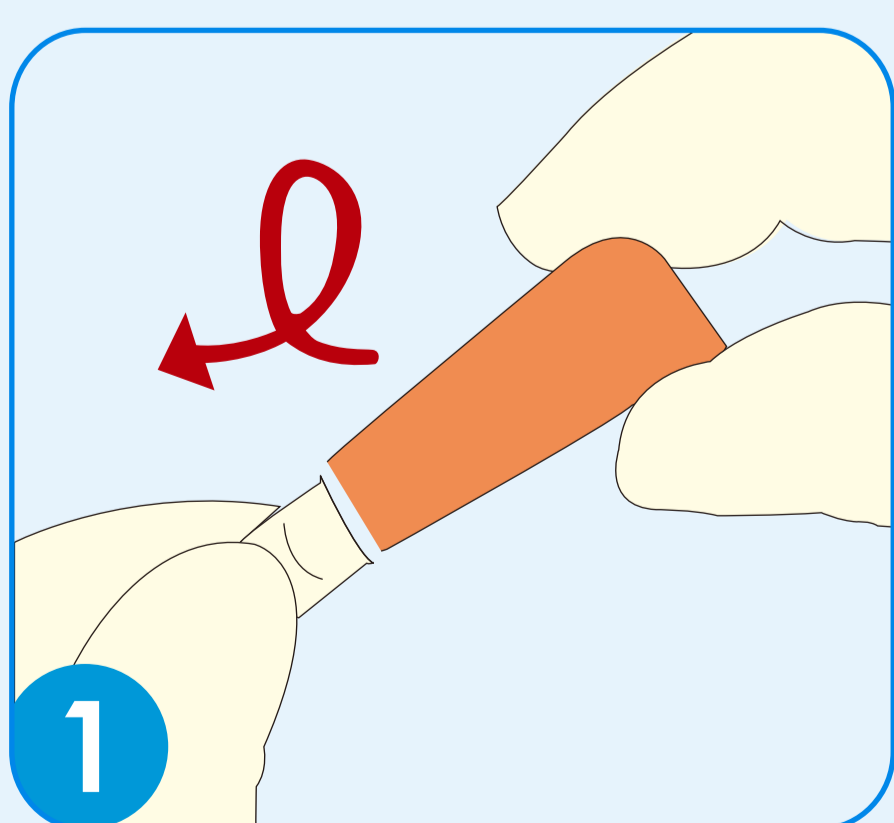

1. Twist and pull the protective cap off. Discard it into the container.

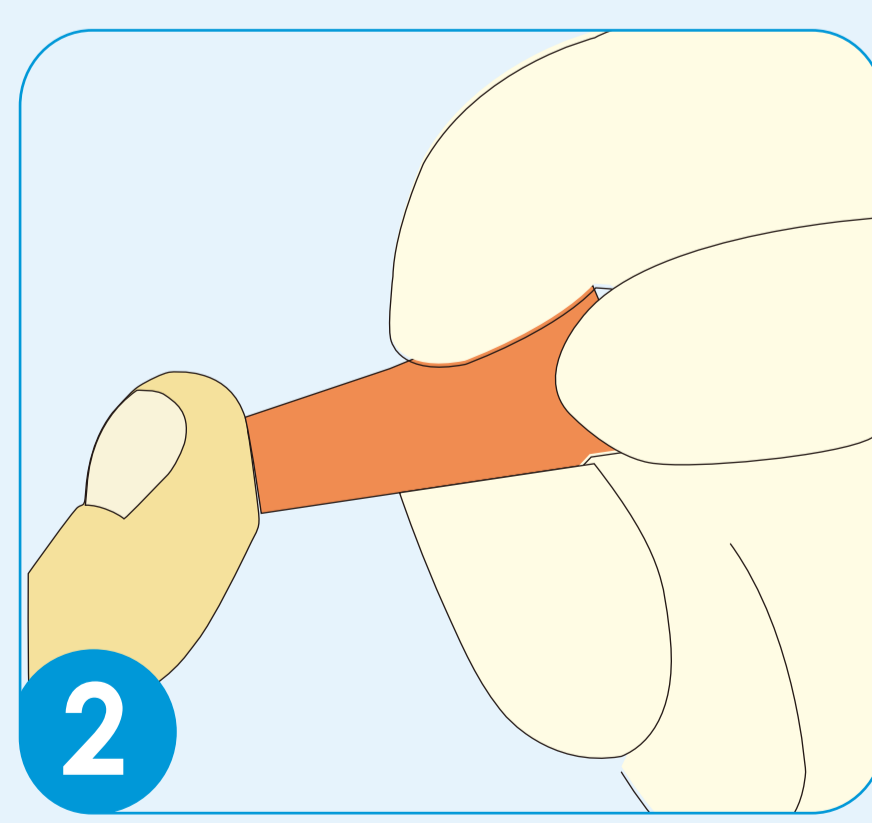

2. Press the lancet against the blood sampling site to get the sample.

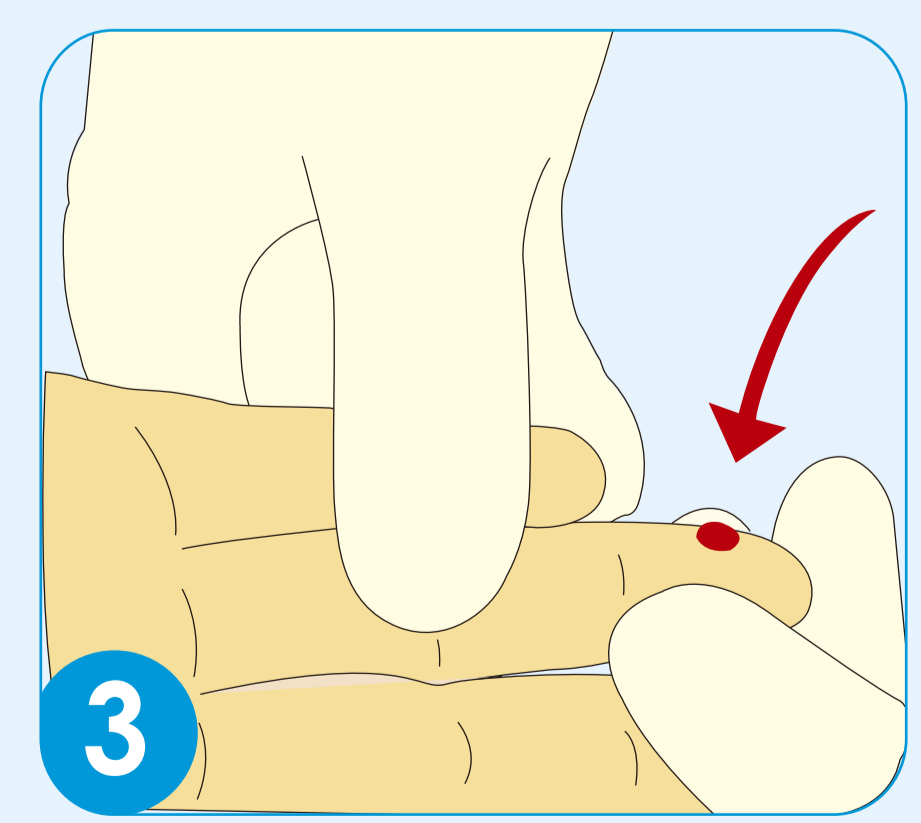

3. Massage the sample site gently to get the sample, and discard the used lancet into a medical waste container.

Note: Please dispose the lancet into appropriate container at once after using.