




*progress in surgery*

## UVDDID HYBRID

Hysteroscopy  
HF SURGERY  
SETS

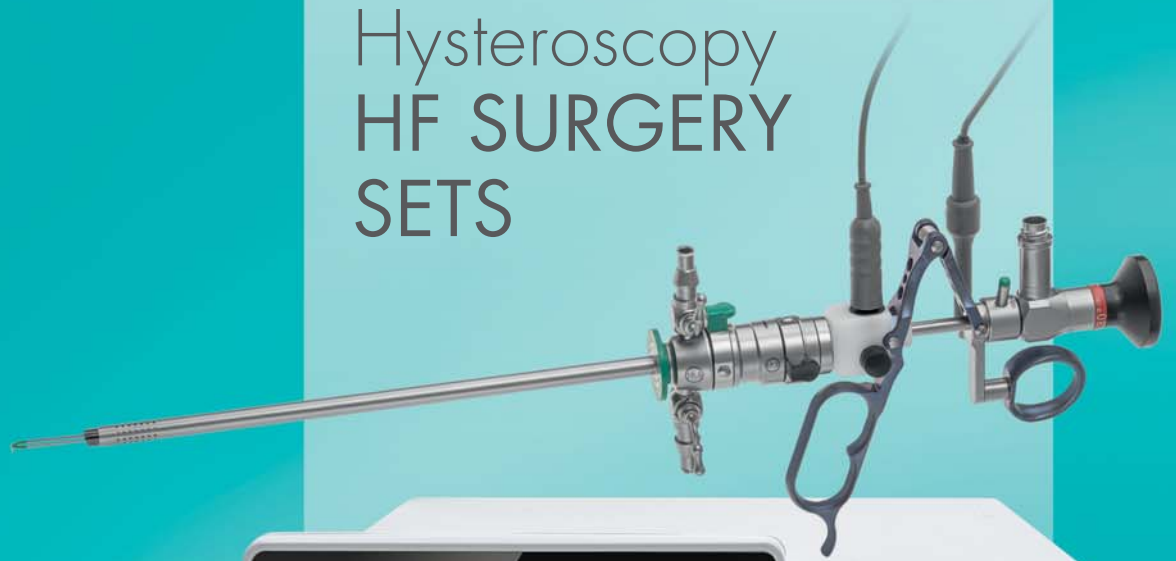

# HYSTEROSCOPY


## PLUG AND PLAY SMART HF SURGERY SETS

The new smart HF surgery combination consisting of G44 unit and endo-hysteroscopic instrumentation provides a perfect variety of plug-and-play modes.

By using of those plug-and-play modes, surgery becomes remarkable easy.

G44 unit is excellent for e.g. bipolar and monopolar resection, vaporization or encucleation. Its increasing safety stops bleeding and reduces operation time.

- Plug and play connectivity makes installation and operation in Hysteroscopy easy with the G44 HF device
- Automatic power regulation
- No configuration is required - but can be made
- Touch display (easy cleaning)
- The OR team will experience state-of-the-art HF surgery with highest simplicity for various Hysteroscopic applications

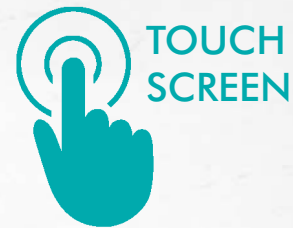

YouTube setup and operation Video on the RZ YouTube Channel:  
<https://youtu.be/Pr-1Er1EO8M>



The smart HYBRID Hysteroscopy Resection Set consists of the HF Generator, the Footswitch, the Instruments, the Neutral Electrode and one bipolar as well as one monopolar cable. Due to the coded SDS-plugs, the surgeon just needs to connect the instrument with the generator and is able to start working without preconfiguring the settings. Herewith we fulfill the state of the art operation standard in the hysteroscopic field.

The RZ HYBRID Resectoscope can be used in the common bipolar as well as in the monopolar way. The Surgeon just needs to change the cable, the rinsing fluid and connect the neutral electrode for monopolar usage. The Resectoscope as well as the Electrodes can be used in an hybrid way. The small 18.5Charr system comes up with numerous advantages like the optimal and painless use with virgo patients.

Even with this small and handy system the surgeon is able to treat all the hysteroscopic pathologies like polyps, myomas, uterine septa, isthmocoele and endometrial ablation. Direct vaginoscopy access is possible and no cervical dilation is needed to proceed time-saving surgeries with 16 Charr. Electrodes and Cold Knife components.

To protect the instruments and extend their lifetime, the instruments are delivered in an optimal fitting storage tray.



G44 Electrosurgical Generator

- Small and versatile generator
- Pre-configured operating modes or free adjustments for special needs while operating
- Simple and intuitive handling and large touchscreen
- With Smart Device Sockets (SDS), G44 recognizes which instrument is connected and uses the appropriate settings
- Changed settings during operation are kept in memory and are recalled for the next procedure
- Automatic neutral electrode monitoring in Monopolar mode - generator stops and signals insufficient connection



Footswitch



Neutral Electrode with cable

18.5 CHARR. **HYBRID** RESECTION SET



351-900-913  
18.5 Charr. HYBRID Resection Set

700-340-400	G44 Electrosurgical Generator 230V*
700-350-405	Footswitch
351-829-030	Hysteroscope, 30°, Ø 2.9mm, autoclavable
351-400-201	Working Element, active
351-000-185	Resectoscope Sheath, rotatable, ceramic tip
351-000-201	Loop Electrode
351-000-203	Knife Electrode
351-000-205	Ball Electrode
351-000-180	Visual Obturator, with 5 Charr. working channel
253-300-999	Sterilisation tubes for electrodes
600-485-555	Stainless Steel Tray, perforated, 480x254x55mm, with lid
600-281-003	Silicone mat, green, 440x230mm
351-400-605	SDS Cable, 4.5m, bipolar
351-501-453	SDS Cable, 4.5m, monopolar
700-503-003	Cable for Neutral Electrode
700-503-000	Disposable Neutral Electrode, pack of 50 pcs.



Complete set in Sterilization Tray

\* also available in 110V: 700-340-410



The smart HYBRID Hysteroscopy Resection Set consists of the HF Generator, the Footswitch, the Instruments, the Neutral Electrode and one bipolar as well as one monopolar cable. Due to the coded SDS-plugs, the surgeon just needs to connect the instrument with the generator and is able to start working without preconfiguring the settings. Herewith we fulfill the state of the art operation standard in the hysteroscopic field.

The RZ HYBRID Resectoscope can be used in the common bipolar as well as in the monopolar way. The Surgeon just needs to change the cable, the rinsing fluid and connect the neutral electrode for monopolar usage. The Resectoscope as well as the Electrodes can be used in an hybrid way.

This 26 Charr. Resectoscope with its classic instrumentation can be used for large-scale myomas, polyps and endometrial ablation.

To protect the instruments and extend their lifetime, the instruments are delivered in an optimal fitting storage tray.



**G44 Electrosurgical Generator**

- Small and versatile generator
- Pre-configured operating modes or free adjustments for special needs while operating
- Simple and intuitive handling and large touchscreen
- With Smart Device Sockets (SDS), G44 recognizes which instrument is connected and uses the appropriate settings
- Changed settings during operation are kept in memory and are recalled for the next procedure
- Automatic neutral electrode monitoring in Monopolar mode - generator stops and signals insufficient connection




Footswitch




Neutral Electrode with cable

26 CHARR. **HYBRID** RESECTION SET

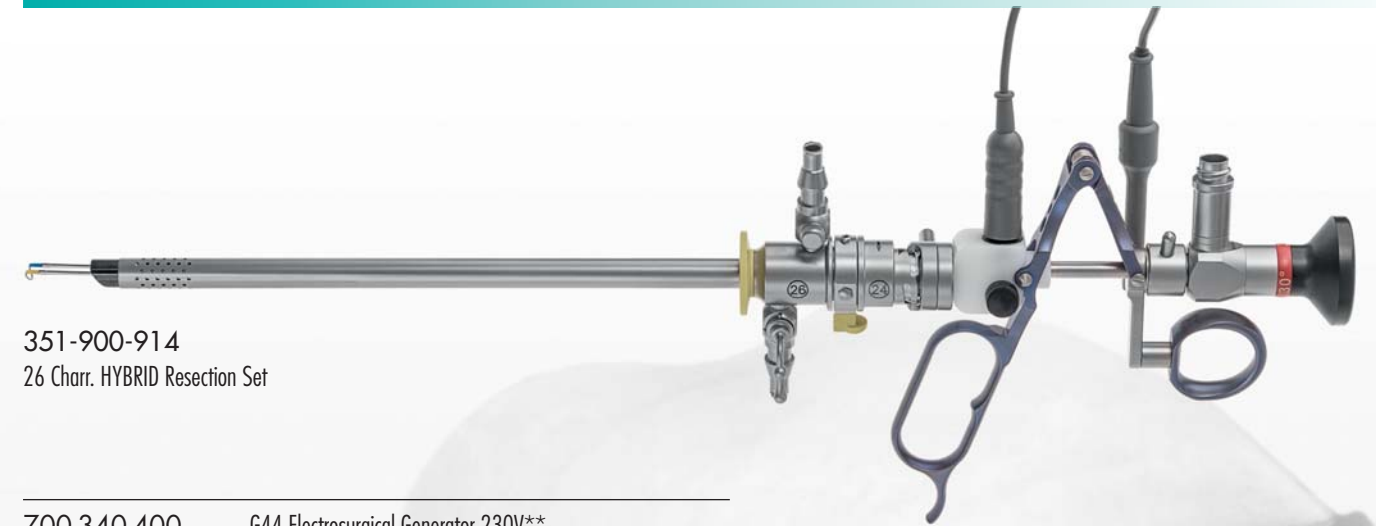

351-900-914  
26 Charr. HYBRID Resection Set

700-340-400	G44 Electrosurgical Generator 230V**
700-350-405	Footswitch
253-804-030	RZ Cystoscope, 30°, Ø 4mm, autoclavable
253-000-401	HYBRID Working Element, active
253-000-256	Resectoscope Sheath, rotatable, ceramic tip
253-400-201	Cutting Electrode
253-400-203	Collins Coagulation Electrode
253-400-205	Ball Electrode
253-000-094	Visual Obturator, with 6 Charr. working channel
253-300-999	Sterilisation tubes for electrodes
600-485-555	Stainless Steel Tray, perforated, 480x254x55mm, with lid
600-281-003	Silicone mat, green, 440x230mm
351-400-605	SDS Cable, 4.5m, bipolar
351-501-453	SDS Cable, 4.5m, monopolar
700-503-003	Cable for Neutral Electrode
700-503-000	Disposable Neutral Electrode, pack of 50 pcs.

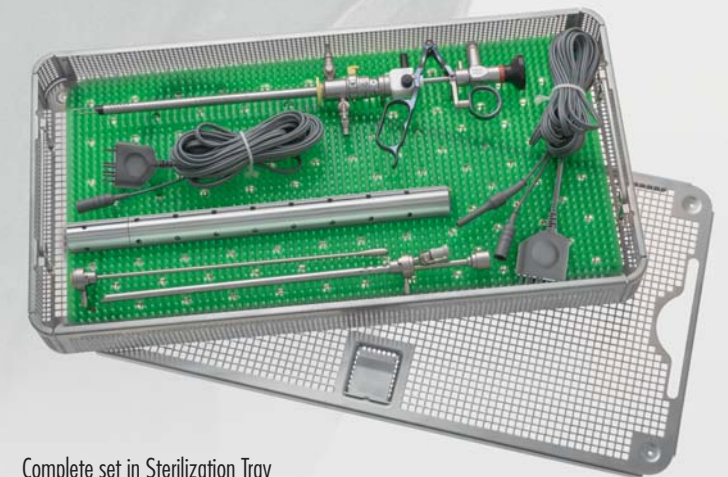

Complete set in Sterilization Tray

\* also available in 110V: 700-340-410



*progress in surgery*

## RZ Medizintechnik GmbH

Unter Hasslen 20  
D-78532 Tuttlingen

Phone: +49 (0)7462/9470-0  
Fax: +49 (0)7462/9470-50

[sales@rz-medizintechnik.com](mailto:sales@rz-medizintechnik.com)  
[www.rz-medizintechnik.com](http://www.rz-medizintechnik.com)

Media Code Number: 021430