

HELMET WITH ZIPPER
FOR CPAP THERAPY
WITH SILICONE COLLAR



HAROL
Your breathing
our mission

FAST-HAR

The helmet **FAST - HAR** is intended to be used for non-invasive ventilation with CPAP therapy in hospital environment. The use of the helmet allows to perform the **therapy in positive pressure with much higher flows** than those used with face masks.

INDICATIONS FOR CLINICAL VENTILATION THERAPY:

Cardiogenic pulmonary oedema - Acute respiratory insufficiency (IRA) - Non-cardiogenic pulmonary oedema (A.L.I. - ARDS) - Chest injury - Atelectasis - Post-operative hypoxemia - Asthma.

» TECHNICAL FEATURES

- The perfect sealing is guaranteed by an ultra-soft **silicone collar** (the only helmet on the market to be equipped with a silicone collar).
- The helmet is equipped with **padded fastening straps** (with neutral Ph) and equipped with a **quick coupling / release device** (with the optional pillow it is possible to use the helmet even without straps).
- A **pocket on the flow inlet connection** directs the flow tangentially to the tent to reduce noise and to avoid direct arrival on the patient's face.
- The **bidirectional anti-suffocation valve** is integrated in the access port and automatically opens when the pressure inside the helmet falls below 2 cmH₂O.
- The transparent tent is made of PVC, the device is LATEX FREE and DEHP FREE.
- Tent volume: 15 lt., device weight: 410 gr.
- The 22F / 22M connections comply with the ISO standard.

» BENEFITS

- EASY TO USE:** the one-piece helmet is practical and easy to put on or take off.
- QUICK ACCESS TO THE PATIENT:** thanks to the convenient zipper, it is possible to carry out face cleaning or bronco-aspiration procedures by a single operator.
- COMFORTABLE:** the optional inflatable pillow further improves patient comfort, provides support to the cervical area, ensures more stability to the helmet, and allows its use without underarm fastening straps.
- PERFECT SEALING:** the exclusive silicone collar ensures the sealing on the neck without having to excessively pull the straps.
- MODULAR SYSTEM:** with this model it is possible to switch from using the helmet with underarm straps to using it without them and vice versa.


We recommend the use of the helmet with devices and equipment suitable for CPAP continuous flow ventilation (minimum flow value 40 l/min.).

THIS IS A DISPOSABLE AND NON-STERILE DEVICE.

Integrated bidirectional anti-suffocation valve

Optional Integrated inflatable pillow

Pressure detection port

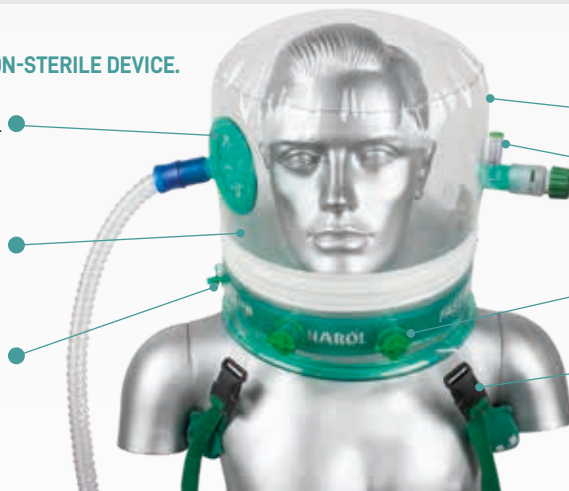

Totally transparent tent

Integrated manometer

Sealed ports for catheters or probes

Padded straps with quick release device

THE HELMET IS AVAILABLE IN VARIOUS SIZES DEPENDING ON THE CIRCUMFERENCE OF THE PATIENT'S NECK.

AVAILABLE PRODUCTS*

ITEM NUMBER	SIZE	NECK CIRCUMFERENCE	PACKAGING
CPAP10451/XS	EXTRA SMALL	28-30 cm	Single packed
CPAP10451/S	SMALL	31-34 cm	Single packed
CPAP10451/M	MEDIUM	35-38 cm	Single packed
CPAP10451/L	LARGE	39-42 cm	Single packed
CPAP10451/XL	EXTRA LARGE	43-45 cm	Single packed
CPAP10451/XXL	2X LARGE	46-49 cm	Single packed
CPAP10451/3XL	3X LARGE	>= 50 cm	Single packed

*PEEP valve and tubing not included and have to be ordered separately. Please refer to the Harol Customer Service for the range of accessories and custom packs.

MADE IN ITALY

HAROL

Via Marcora 53, 20097
San Donato Milanese (MI) - ITALY
+39 02 5272788

info@harol.it
www.harol.it

